

Automatic filling and closing machine Minican pharma Equipment:





Indexing table with 8, 12 or 16 stations

1-4 Filling stations

Application station for rubber stoppers or droppers

Application station for alu caps or screw caps

Crimping station or screw capping station

Optionals:

Laminar flow

Check weigher

Labelling machine

Ink Jet Coder

Bottle washing machine

and others

bottles and/or vials

Tel.: **43/2235/877 22-0

Fax: **43/2235/877 22-16



Dosing systems:

mass flow meter

magnetic inductive flow meter

turbine flow meter

peristaltic dosing pump

piston dosing pump

and others



Filling ranges: 0,1 – 2000 ml (optional up to 3000 ml)

Tel.: **43/2235/877 22-0

Fax: **43/2235/877 22-16

Output ranges: 20 40, 60, 80 and 120 bottles/min

Product contacting parts: AISI 316 L, others to be agreed

Optionals:

Labelling units

Check weigher

Carton packers

Complete data processing system

Inspection systems

Validation packages

Engineering support for turnkey projects







The massflow-meter

one good option for the dosing system



Filling hose

← Flow speed regulation

← Flow meter

Tel.: **43/2235/877 22-0

Fax: **43/2235/877 22-16

Product infeed (0,5-2 bar)

Filling nozzle with mounted CIP/SIP waste hose



Filling nozzle (reaching inside CIP/SIP fitting)

CIP/SIP fitting with temperature sensor for streaming temperature

Waste hose for cleaning medium and steam



The peristaltic pump another brilliant option for the dosing system



Quick-release filling needles



Tel.: **43/2235/877 22-0

Fax: **43/2235/877 22-16





Application with rubber stopper/crimp cap and screw cap




Application with labeller

Tel.: **43/2235/877 22-0 Fax: **43/2235/877 22-16



Example of interface settings

(in reference to the machine handbook)

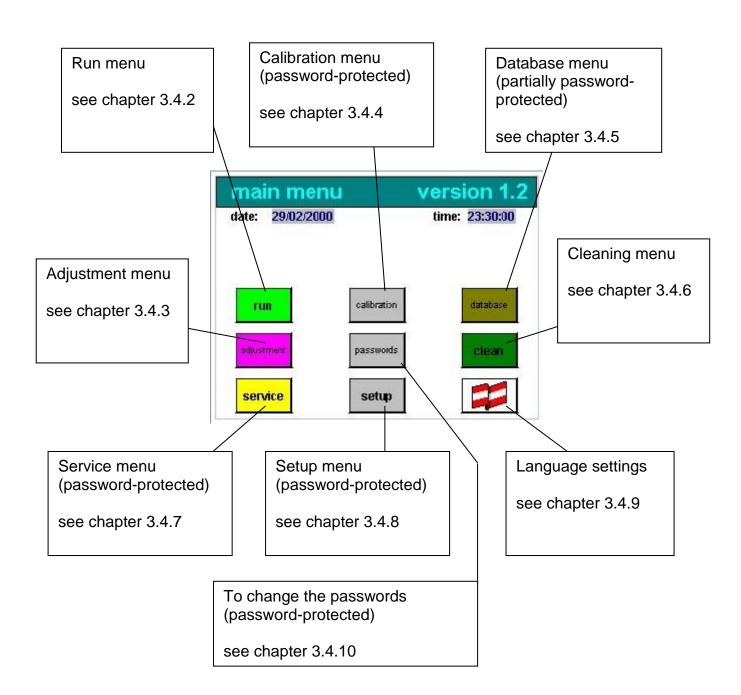

Tel.: **43/2235/877 22-0

Fax: **43/2235/877 22-16